

Infrared Chemical Imaging: Semi-Quantitative Analysis of Components

Jonathan Tarr, Koichi Nishikida, Federico Izzia
 Thermo Electron Corporation
 jonathan.tarr@thermo.com

With new image analysis software for infrared spectroscopy, semi-quantitative studies can be performed without the use of standards providing rapid measurements of component concentrations when an approximate value is suitable or when quantitative standards are not available. Infrared spectroscopy has been used for many years to identify the composition of organic and other materials. Single or multi-component samples can be analyzed and their chemical structures determined by

comparing their absorbance peaks to known spectral libraries or fundamental molecular transitions. This technique is used throughout the scientific community and provides rich information.

More recently, chemical mapping or imaging has been developed to measure both the compositions and spatial distributions of components in a sample. This gives two-dimensional information of compositions and provides more useful data when relative amounts or positions are important. In order to perform quantitative analysis, standards must be analyzed and the instrument's response characterized to create a calibration curve for each component of interest. This work is tedious as the analyst must spend time preparing the standards, analyzing them in addition to creating the mathematical equations that will be used to model the unknown samples of interest.

In the following studies, semi-quantitative results are obtained using chemical images and powerful image analysis software. The images are first binarized and then the relative areas of components are calculated to give a percent of the total area. This number is in direct correlation to the actual quantitative values, yet no standards are needed. A more accurate value can be obtained by taking into account the densities of the components, yielding a weight percent instead of an area percent concentration. For many purposes, the semi-quantitative result is good enough and is close to the actual value within a few percent. However, in some analyses, quantitative standards are difficult or impossible to obtain. Two case studies are presented: A generic pharmaceutical tablet and a micro-tomed section of a coffee bean.

Generic Pain-Relief Tablet

A pharmaceutical tablet containing acetaminophen and binder material was analyzed using a Thermo Electron Corporation Nicolet™ Continuum™ XL FT-IR imaging microscope. The data was collected in reflectance mode, measuring the absorbance of reflected infrared energy. The tablet was approximately 8x15 mm in size. A representative map measuring approximately 1.5x1.5 mm was collected in the spectral range of 4000-800 cm⁻¹. The chemical images, video image and representative spectra of acetaminophen and the cellulosic binder material are shown in figure 1. The chemical images, or maps, correlate to the spectral intensities of the data at 1655 and 1066 cm⁻¹ for acetaminophen and cellulose, respectively. Binarized images were created using OMNIC™ Atlas™ image analysis software. The areas in square microns of the two components were calculated by the software. The relative percentages of

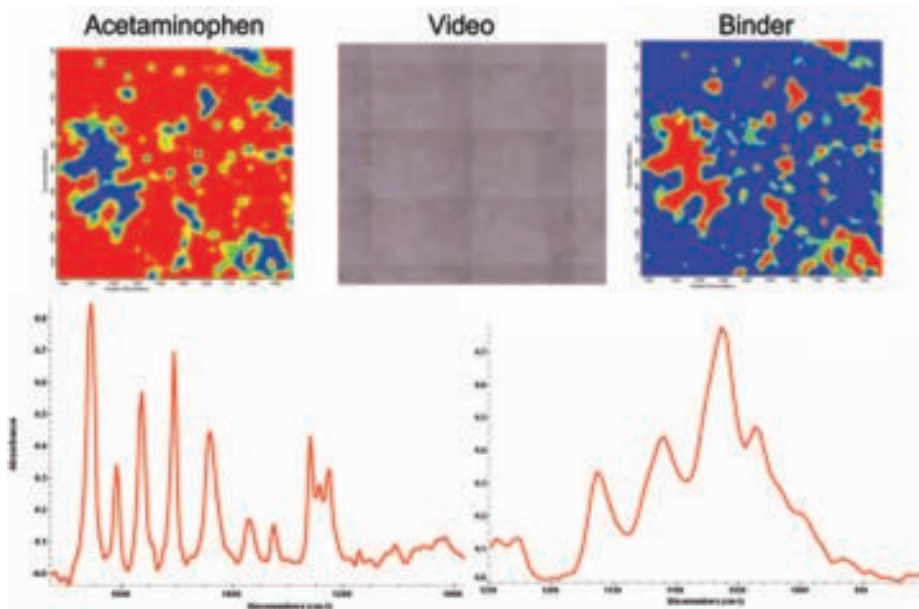


Figure 1. Video image, representative component spectra and corresponding chemical maps of generic pain-relief tablet.


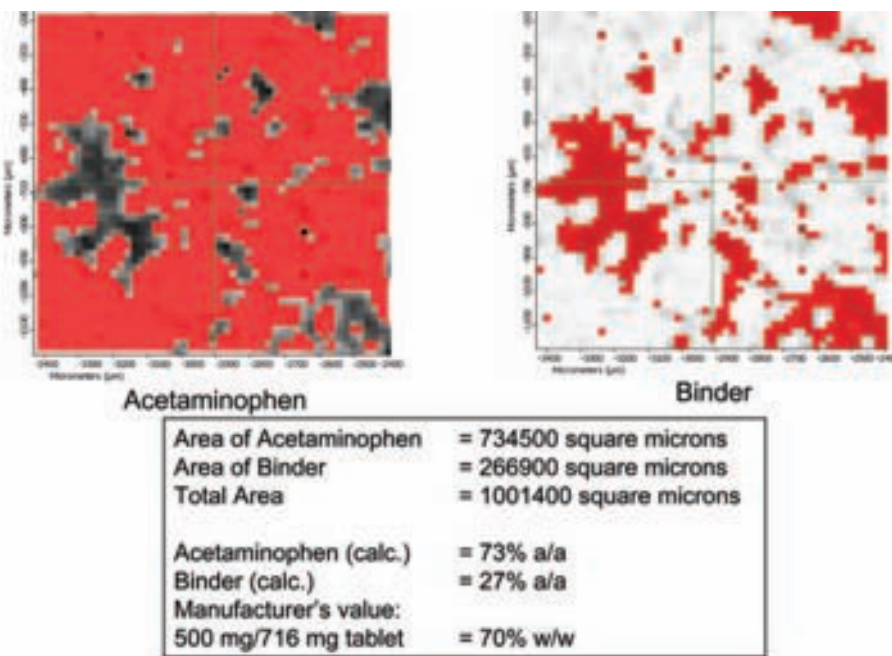


Figure 2. Image analysis of pain-relief tablet showing the binarized chemical maps and quantitative results table.

The Next Step... Innovative Solutions for Nano-Technology

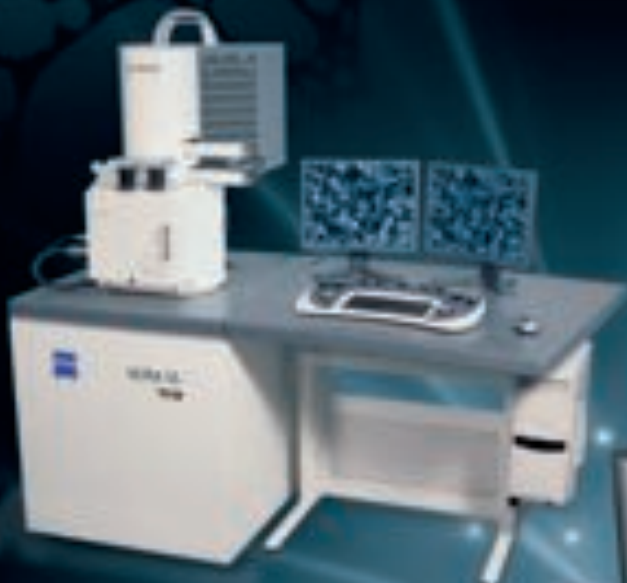

ULTRA FESEM

An extension of our leading SUPRA™ series with new In-column EsB detector for superb high resolution imaging and unmatched compositional data.



EVO® SEM

Tungsten or LaB₆ system with superior imaging quality, class leading analytical capabilities and extened pressure option.



LIBRA® EFTEM

Unique In-column OMEGA filter combined with Koehler illumination for unrivalled flexibility in imaging and analysis.

Enabling the Nano-Age World®

Carl Zeiss SMT Inc
Nano Technology Systems Division
One Zeiss Drive
Thornwood, New York 10594
USA

Tel. +1 914/747 7700
Fax +1 914/681 7443
info-usa@smt.zeiss.com
www.smt.zeiss.com/nts



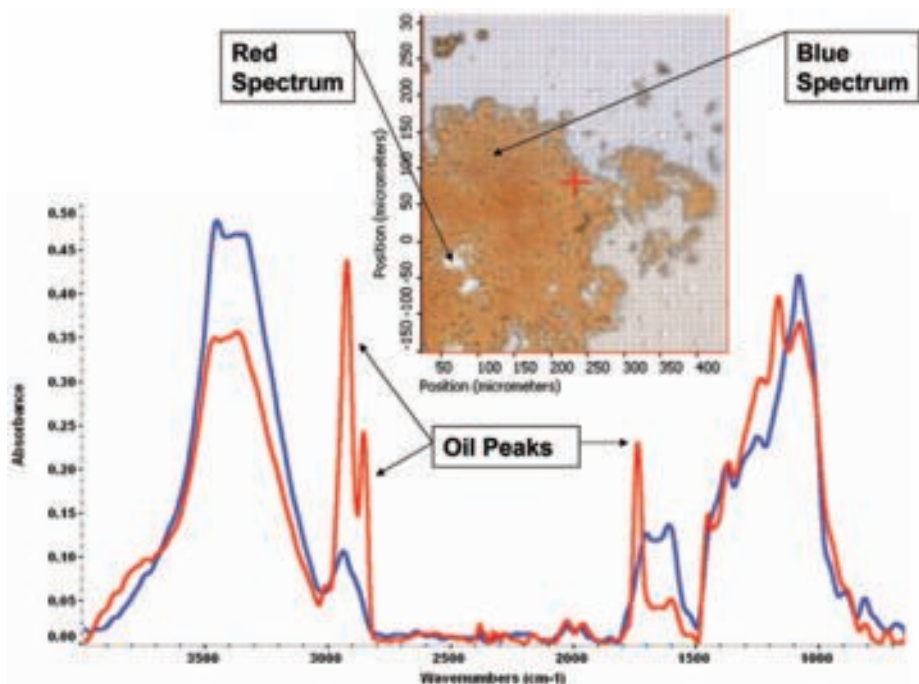


Figure 3. Video image of coffee bean section and spectra from the two points referenced. The blue spectrum shows a cellulose structure and the red spectrum shows contributions from cellulose and oils.

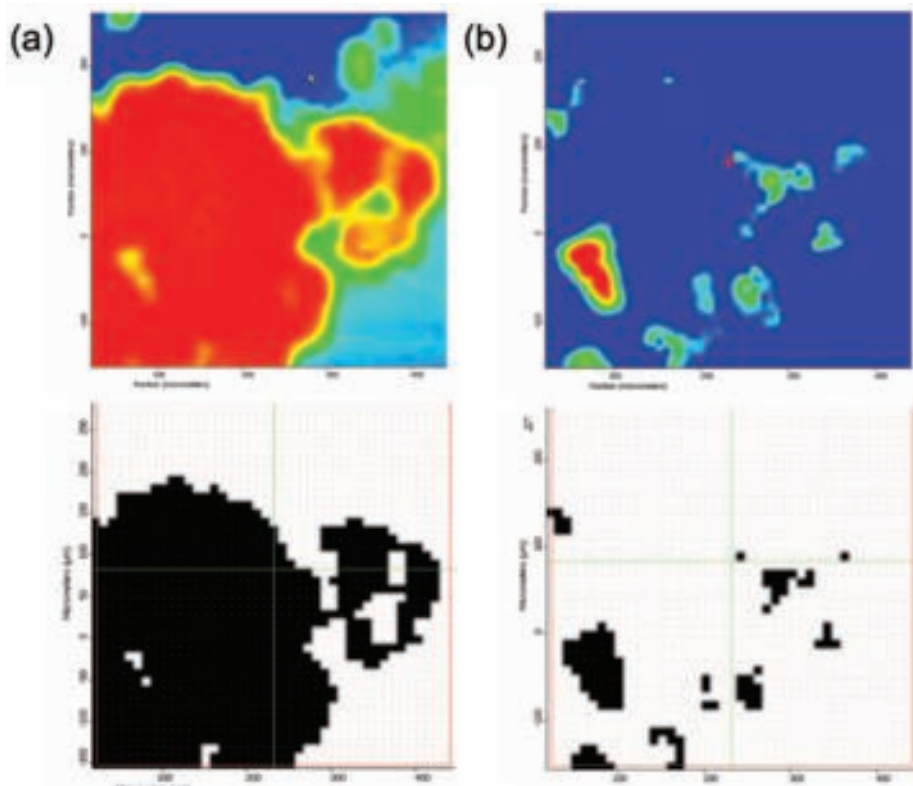


Figure 4. Spectral correlation maps and corresponding binarized data for cellulose (a) and oil-containing (b) areas.

Area of Cellulose	= 95911 square microns
Area of Oils	= 10592 square microns
Total Area	= 106503 square microns
Cellulose	= 90% a/a
Oils	= 10% a/a

Figure 5. Quantitative results of image analysis applied to the binarized coffee bean data.

the total area were then calculated for each component. The tablet was weighed using an analytical balance. The manufacturer listed the amount of acetaminophen in the tablet and the rest of the mass was assumed to be binder material. Using this information, a theoretical weight percent value was determined for the active ingredient and binder. These values were compared to those obtained from the image analysis software. The binarized images and results are given in figure 2. The experimental values obtained are within 3% of the theoretical values.

Coffee Bean Section

A fresh roasted coffee bean was micro-tomed into a 20-micron section. The section was placed in a diamond compression cell and analyzed using infrared transmission spectroscopy. The video image and spectra from two points, A and B, on the map are presented in figure 3. The spectrum at point A shows a typical cellulosic structure and the spectrum at point B (the large clear area in the video image) shows both cellulose and vegetable oil components. Spectral correlation maps of the cellulose spectrum and the oily spectrum as well as their respective binarized images are given in figure 4. The two correlation maps were obtained by comparing the two spectra, A and B, to all of the spectra in the data set and mathematically assigning a correlation coefficient to each. The binarized data was analyzed and the total area of each component was calculated as a percentage of the total area. Figure 5 gives the results of the analysis. The area of oils was determined to be 10% of the total area.

Conclusions

Infrared imaging in combination with image analysis allows the composition of heterogeneous materials to be determined without the need for external standards or calibration. A high-contrast chemical image can be created based on the IR spectral signature that is unique to each component. Area analysis of the image then yields the quantitative composition of the sample. In other words, infrared imaging provides information regarding the identity, quantity and location of components within complex materials. In the case of the pharmaceutical tablet image analysis provided a close estimate of the known concentrations of components. With the coffee bean section, there was no known compositional data for reference. Image analysis provided an estimate of the oil concentration that can be used to predict the flavor content of the coffee. ■

NEW **AUTOPROBE™ 300**

**HIGH THROUGHPUT HIGH SUCCESS RATE
SAVES FIB TIME**

The **NEW** AutoProbe™ 300 is the backbone of any high throughput TEM sample preparation strategy.

• AutoProbe™

8th generation nanomanipulator based on Omniprobe's patented technology

Automated, in-situ, multi-purpose system capable of :

- TEM sample lift-out
- Electrical testing
- Nano-mechanical testing
- Charge neutralization

• IN-SITU PROBE TIP EXCHANGE

Hardware and software package for in-situ probe tip exchange

- Custom probe shaft
- Low cost consumable probe tips
- Multi-tip cassette for tip transfer from FIB to Short-Cut™
- Contact detection with built-in strain gauge

• FAST AND EASY TEM SAMPLE PREP

Short-Cut™ tool enables ex-situ sample transfer from probe tip to 3mm TEM grid

- Bypasses the in-situ grid attach step
- Releases FIB bottleneck and optimizes FIB utilization
- Low cost consumable TEM grid

Extensive Patent Portfolio

CE and UL certified

Omniprobe, Inc.
10410 Miller Road
Dallas, TX 75238
Phone: 214.572.6800
Fax: 214.572.6801

omniprobe
Extending FIB Capabilities
www.omniprobe.com

For Sales: sales@omniprobe.com
For Support: customer-service@omniprobe.com

European Support:
SEM-FIB Solutions
+44 1993 771203
tony.kendrick@sem-fib.com
Oxford, United Kingdom