

## Progress Towards More Realistic In-Situ Microscopy Observations

A. Howie

Cavendish Laboratory, University of Cambridge, UK

As progress indicators in electron microscopy, advances in spatial resolution and in spectroscopy probably attract most frequent attention. Improved user-friendliness has also been significant even when judged in comparison with scanned probe microscopy. Evidence for developments in in-situ microscopy at least equally impressive can be found by comparing the relevant sections of the book by Hirsch et al. [1] and a more recent compilation [2].

Peter Hirsch's research group swiftly discovered the power and frustrations of in-situ microscopy. The cine film of dislocation motion observed in the earliest diffraction contrast studies [3,4], and attributed to the thermal stresses generated by the electron beam, was extremely effective in convincing the wider community that it was indeed dislocations that were being observed. In a remarkably successful in-situ hot stage experiment applied to a simple problem, the shrinkage of dislocation loops due to migration of vacancies to the foil surfaces was followed and fitted in detail to a quantitative model [5]. Any consequent feelings of euphoria were probably moderated when, in trying to employ the same equipment more ambitiously to study alloy phase transitions, the dominating effect of diffusion and reactions at the foil surfaces was observed [6]. Another hot stage project that was tackled with mixed success was to follow the recrystallisation process of grain boundary motion in heavily worked silver foils [7]. Successful cine photography was defeated by the tendency for recrystallisation to take place abruptly outside the field of view. Observations of dislocation motion under more controlled conditions of plastic deformation as well as in fracture have continued but the focus below is on surface science and catalysis work.

Despite the poor vacuum conditions, the power of TEM in studies of metal epitaxy was apparent [8] and encouraged Hirsch and Pashley to enlist the expertise of Ugo Valdre in the construction of a separately pumped UHV stage unit [9]. This system was used for several years in studies of both island and pseudomorphic growth as well as in metal oxidation [10] where much of the work has been confirmed in more modern equipment [11] capable of revealing the surface reconstructions associated with the oxide nuclei. For clean surface studies, these TEM methods are increasingly challenged by other techniques, including reflection electron microscopy (which first revealed the nucleation of surface reconstructions at surface steps [12]), photoemission electron microscopy (PEEM)[13] and, most particularly, scanned probe microscopy (SPM), whose popularity with the wider surface science community vividly illustrates the importance of user-friendliness. These other approaches have been compared [14] and offer varying levels of spatial resolution, nanoanalysis and imaging speed. They generally employ bulk sample surfaces rather than thin film surfaces with advantages for specimen preparation and surface cleaning but, in some cases, with disadvantages for ionisation damage and beam charging.

The drive towards in-situ microscopy of catalytic processes received early stimulation from some spectacular TEM observations

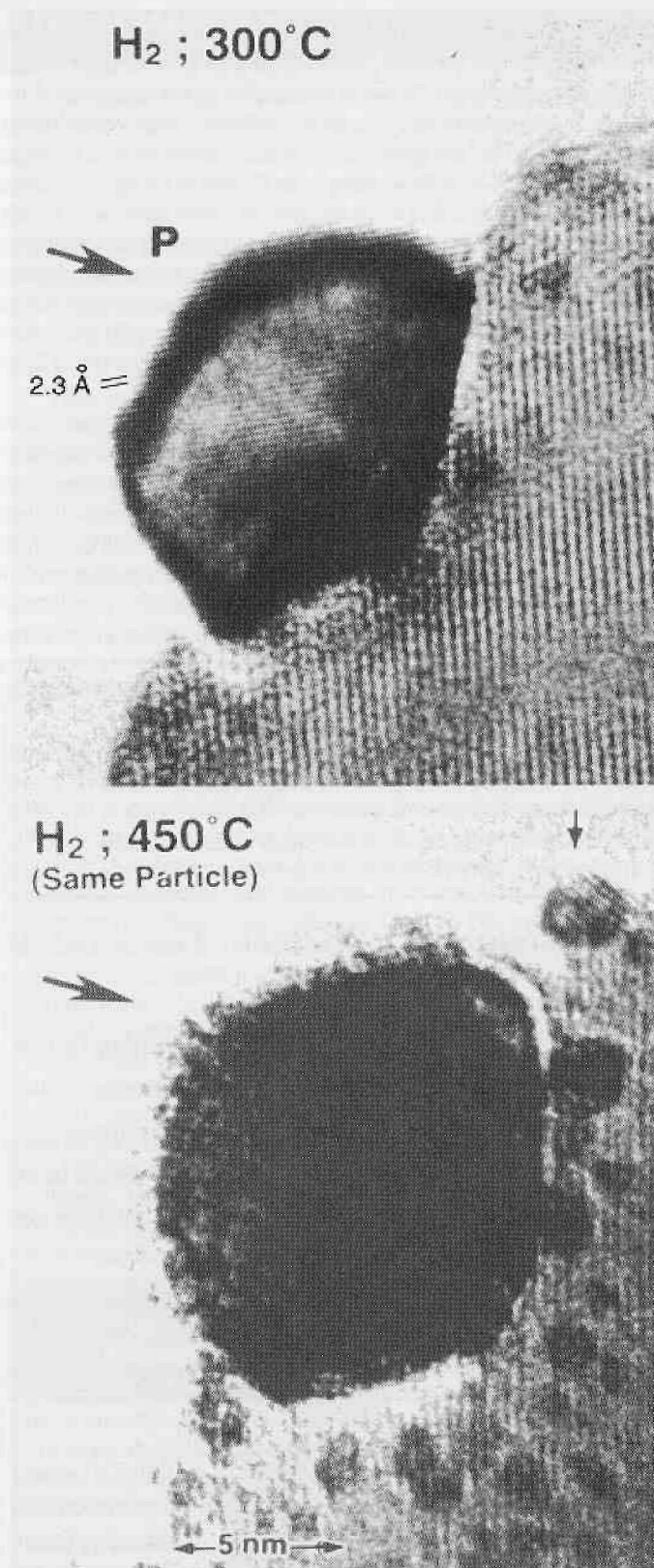


Fig.1. Pt/titania: In situ ETEM at the atomic level, showing strong metal-support interactions: Top: Pt particle at 300°C in H<sub>2</sub>. Bottom: The same particle at 450°C in H<sub>2</sub>. The image shows an overlayer of titania on the particle and the formation of small clusters with ~0.2 nm lattice spacings (arrowed). Taken from ref. [19] courtesy P.L. Gai and Kluwer Academic/Plenum Publishers.



[15] but places extreme demands if realistic pressures and temperatures are to be achieved. Peter Hirsch's long lasting influence here can be traced to the Oxford project for high voltage in-situ microscopy [16]. In some sense this is the ancestor of the much more ambitious medium voltage, high resolution, environmental TEM equipment designed and developed by E. Boyes and P.L. Gai [17] after their move from Oxford to Dupont and used with remarkable success in many catalyst investigations [2, 18, 19]. Fig. 1 shows fine detail in an example of the strong metal-support interaction. A later version of this microscope, installed at Haldor Topsoe A/S, recently yielded intriguing results on the distribution of Ba promoter in Ru ammonia synthesis catalysts [20] and other remarkable in-situ TEM studies of Ziegler-Natta catalysts have been reported with related equipment [21].

Environmental SEM and PEEM both offer some competition to TEM methods for work under gas pressure but cannot compete in terms of spatial resolution. The most serious potential competitor here is undoubtedly SPM but for some reason the challenge has been slow to develop. Recent STM work in the few bar range [22, 23] could surely be extended to higher pressures or even liquid operation conditions on catalysts (provided any oxide support is in the form of a thin coating on metal [24]). Atomic force microscopy and some other forms of SPM using optical methods, such as sum frequency generation [25] could also be potentially powerful methods for imaging catalytic reactions under gas pressure.

A general and incompletely resolved problem for all in-situ work, other than in SPM, is the possible influence of the electron beam. Atomic displacements are often manifest, for instance in HREM profile imaging of surfaces, but in other cases may still affect the course of any dynamic process even when they are not directly

visible. In this respect we have not greatly advanced beyond the early observations of beam-induced dislocation motion [26]. ■

## References

- [1] P.B. Hirsch et al., *Electron Microscopy of Thin Crystals*, Butterworth, London, 1965.
- [2] P.L. Gai (ed), *In-situ Microscopy in Materials Research*, Kluwer, Boston, 1997.
- [3] P.B. Hirsch et al., *Phil. Mag.* 1 (1956) 677.
- [4] M.J. Whelan et al., *Proc. Roy. Soc. A* 240 (1957) 524.
- [5] J. Silcox and M.J. Whelan, *Phil. Mag.* 5 (1960) 1.
- [6] G. Thomas and M.J. Whelan, *Phil. Mag.* 6 (1961) 1103.
- [7] J.E. Bailey, *Phil. Mag.* 5 (1960) 833.
- [8] D.W. Pashley et al., *Phil. Mag.* 10 (1964) 127.
- [9] U. Valdre et al., *J. Phys. E* 3 (1970) 501.
- [10] D.A. Goulden, *Phil. Mag.* 33 (1976) 393.
- [11] J.C. Yang et al., *Microscopy Microanal.* 7 (2001) 486.
- [12] K. Yagi, *Surf. Sci. Rep.* 17 (1993) 305.
- [13] E. Bauer, *J. Electron Spectrosc. & Rel. Phenom.* 114 (2001) 975.
- [14] E.D. Williams and L.D. Marks, *Crit. Rev. in Surf. Chem.* 5 (1995) 275.
- [15] R.T.K. Baker and R.J. Waite, *J. Catal.* 37 (1975) 101.
- [16] P.L. Gai and M.J. Goringe, *Proc 39th EMSA Conf.*, San Francisco Press, 1981, 68.
- [17] E.D. Boyes and P.L. Gai, *Ultramicrosc.* 67 (1997) 219.
- [18] P.L. Gai, *Advanced Mater.* 10 (1998) 1259.
- [19] P.L. Gai, *Topics in Catalysis* 8 (1999) 97.
- [20] T.W. Hansen et al., *Science* 294 (2001) 1508.
- [21] V.P. Oleshko, P. Crozier, R.D. Cantrell and A.D. Westwood, *J. Electr. Microsc.* 51 (2002) S27.
- [22] B.L. Weeks et al., *Rev. Sci. Instrum.* 71 (2000) 3777.
- [23] L. Osterland et al., *Phys. Rev. Lett.* 86 (2001) 460.
- [24] H. Over et al., *Science* 287 (2000) 1474.
- [25] T. Dellwig et al., *Phys. Rev. Lett.* 85 (2000) 776.
- [26] This article is slightly expanded from one prepared for the Quebec MSA meeting's in-situ microscopy symposium honoring Sir Peter Hirsch. I am very grateful to Dr P.L. Gai, one of the organizers, for much helpful advice and for making available the figures used here.

## Introducing EDS2000i, The Most Advanced Microanalysis System Today!

- 100% Integrated inside the SEM operating system
- Apply the EDS tools on the live SEM Image
- Analyze samples, faster, with greater accuracy
- Save all EDS data along with the SEM image in one file
- Integrated Stage control for advanced particle analysis



IXRF's Integrated EDS systems combine the SEM and IXRF EDS user interfaces into a single application. The EDS Toolbar is accessed directly from the SEM user interface allowing advanced microanalysis features to be performed directly on the live SEM image. Spectra, X-ray Maps, and Linescans can be directly acquired from the selected regions of the SEM image with a single *click* on the EDS Toolbar. The EDS and SEM user interfaces combine to provide one seamless microanalysis tool.



IXRF Systems, Inc., Houston, TX, USA, [www.ixrfsystems.com](http://www.ixrfsystems.com)  
Ph:(281) 286-6485 Fax:(281)286-2660



## COMING EVENTS


- ✓ **Microscopy and Microanalysis 2002**  
August 4-8, 2002, Quebec City, Quebec, Canada  
[www.msa.microscopy.com](http://www.msa.microscopy.com)
- ✓ **7th Intl. Conf. on Near Field Optics & Related Techs.**  
August 11-15, 2002, Rochester, NY, USA  
[www.optics.rochester.edu](http://www.optics.rochester.edu)
- ✓ **Microscopical Society of Ireland's Annual Symp.**  
August 28-30, 2002, Galway, Ireland  
[www.nuigalway.ie](http://www.nuigalway.ie)
- ✓ **XVth Intl. Congress on Electron Microscopy**  
September 1-6, 2002, Durban, SA  
[www.ru.ac.za/affiliates/emu/icem2002.htm](http://www.ru.ac.za/affiliates/emu/icem2002.htm)
- ✓ **BIOMED 2002**  
September 19-22, 2002, Kemer Antalya, Turkey  
[www.biomed2002.ankara.edu.tr](http://www.biomed2002.ankara.edu.tr)
- ✓ **National Society for Histotechnology**  
September 28 - October 3, 2002, Long Beach, CA  
[histo@nsh.org](mailto:histo@nsh.org)
- ✓ **American Society of Human Genetics**  
October 15-19, 2002, Baltimore, MD  
1 (301) 571-1825
- ✓ **American Society for Cell Biology**  
December 14-18, 2002, San Francisco, CA  
[www.acsb.org](http://www.acsb.org)
- ✓ **Microscopy NZ Conference**  
February 10-14, 2003, Wellington, NZ  
[peter.smith@agresearch.co.nz](mailto:peter.smith@agresearch.co.nz)
- ✓ **PITTCON 2003**  
March 9-14, 2003, Orlando, FL  
[www.pittcon.org](http://www.pittcon.org)
- ✓ **Scanning 2003**  
May 3-5, 2003, San Diego, CA  
[www.scanning.org](http://www.scanning.org)
- ✓ **Microscopical Society of Canada**  
June 4-6, 2003, Vancouver, BC, Canada  
[ech@unixg.ubc.ca](mailto:ech@unixg.ubc.ca)
- ✓ **Microscopy and Microanalysis 2003**  
August 3-7, 2003, San Antonio, TX  
[www.msa.microscopy.com](http://www.msa.microscopy.com)
- ✓ **American Society for Cell Biology 2003**  
December 13-17, 2003, San Francisco, CA  
[www.ascb.org](http://www.ascb.org)
- ✓ **Microscopy and Microanalysis 2004**  
August 1-5, 2004, Savannah, GA  
[www.msa.microscopy.com](http://www.msa.microscopy.com)
- ✓ **Microscopy and Microanalysis 2005**  
July 31- August 4, 2005, Honolulu, HA  
[www.msa.microscopy.com](http://www.msa.microscopy.com)
- ✓ **Microscopy and Microanalysis 2006**  
August 6-10, 2006, Chicago, IL  
[www.msa.microscopy.com](http://www.msa.microscopy.com)

Please check the "Calendar of Meetings and Courses" in the MSA journal "Microscopy and Microanalysis" for more details

# MRS-4

*A ISO-9000 and ISO Guide-25 Standard for Microscopy  
Calibrate from 10X to 200,000X*

This is our third generation, traceable, magnification reference standard for all types of microscopy (SEM, Optical, STM, AFM, etc.). The MRS-4 has multiple X and Y pitch patterns that range from  $\frac{1}{2} \mu\text{m}$  ( $\pm 0.045 \mu\text{m}$ ) to  $500 \mu\text{m}$  ( $\pm 0.1 \mu\text{m}$ ) and a 6 mm ruler with 1  $\mu\text{m}$  increments.

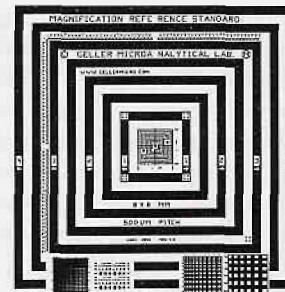

Visit our website and send for our free resources guide.



**GELLER**  
MICROANALYTICAL  
LABORATORY

426e BOSTON STREET (RT. 1) \* TOPSFIELD, MA 01983  
978.887-7000 \* 978.887-6671 \* [jg@gellermicro.com](mailto:jg@gellermicro.com)  
<http://www.gellermicro.com>



## EDAX Clarification

The May/June issue of Microscopy Today published an article by Dr. David Dingley, "Fast Orientation Imaging Microscopy". We wish to clarify that Dr. Dingley is Chief Scientist at TSL; the article was based on work performed at TSL using the company's Orientation Imaging Microscopy, or OIM, software. Dr. Dingley pioneered automated Electron Backscatter Diffraction (EBSD) in the SEM and was a cofounder of TSL, which is now a subsidiary of EDAX.

Orientation Imaging Microscopy uses EBSD to map the spatial distribution of crystallographic orientation and quantify information about a sample's microstructural characteristics. The technique has been extended, under the product name Delphi, to allow identification of unknown crystallographic phases, similar to X-ray diffraction. Capabilities similar to OIM have also been adapted for use in the TEM. This new development is the first fully Automated Crystallography system for TEM (ACT). To learn more, please visit our website at [www.tsl-oim.com](http://www.tsl-oim.com) or e-mail us at [tsl.info@ametec.com](mailto:tsl.info@ametec.com)