

**Helios Software, Inc.**

Charles E. Rose  
P.O. Box 44-1215  
W. Somerville, MA 02144

Tel. 617-776-2943  
FAX 802-785-2042

- MacPDF™ is a Macintosh Search/Retrieval program for the PDF-2. It is a product of Helios Software and is now in use in Japan, Europe and the U.S.A.
- MacPDF now offers software localization for foreign language environments.
- Full Boolean capabilities available December 1992.


*Powder Diffraction* makes this page available for brief announcements of services offered by its readers or by other interested persons.

- Space Cost \$150 (no agency discount). This space will permit the display of a business card or other similar size text. The business card will be reproduced at no charge.
- Free space: Up to 12 lines of text immediately below the announcement.

For details write or call Mary M. Rossi, Advertising Manager, *Powder Diffraction*. ICDD, 1601 Park Lane, Swarthmore, Pennsylvania 19081  
215-328-9403 FAX 215-328-2503

Company	Page Number
Charles Supper Company	XIV
Charles Supper Company	247
Dapple Systems	XII
Fein-Marquart Associates, Inc.	I
Gem Dugout	247
ICDD – British Crystallographic Assoc.	X
ICDD – Clinics	X
ICDD – For 1992/1993	XVII
ICDD – Metals and Alloys	XIV
ICDD – PC-PDF CD-ROM	XV
IC Laboratories, Inc.	XIV
inel Instrumentation Electronique	IV
Ital Structures	XII
KeveX Instruments	XI
Materials Data, Inc.	IX
Omni Instruments, Inc.	VI
Omni Instruments, Inc.	VII
Philips Electronic Instruments, Inc.	Back Cover
Rigaku, U.S.A.	Inside Back Cover
Scintag, Inc.	V
Seifert X-Ray Corporation	XVI
Siemens	Inside Front Cover
Sietronics Pty. Ltd.	XIII
SPEX Industries, Inc.	II

For 1992/93 from ICDD...

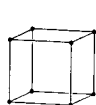

**Set 42 PDF**

In all current media

CD-ROM & Mag Tape... Microfiche ... Book


Over 2000 new &/or updated entries in PDF-1 and PDF-2†



**EDD Electron Diffraction Data Base**

(NIST/Sandia/ICDD)

Crystallographic & chemical data on over 71,000 crystalline materials.†

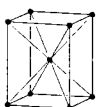

**NIST Crystal Data File**


Crystallographic & chemical data on over 170,000 crystalline materials\*


†Available for IBM-PC compatibles, VAX & Macintosh

\*Available for IBM-PC compatibles & VAX

International Centre for Diffraction Data

1601 Park Lane, Swarthmore, PA 19081 USA (215) 328-9400 Telex: 847170

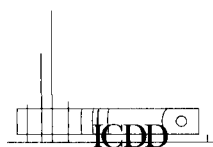

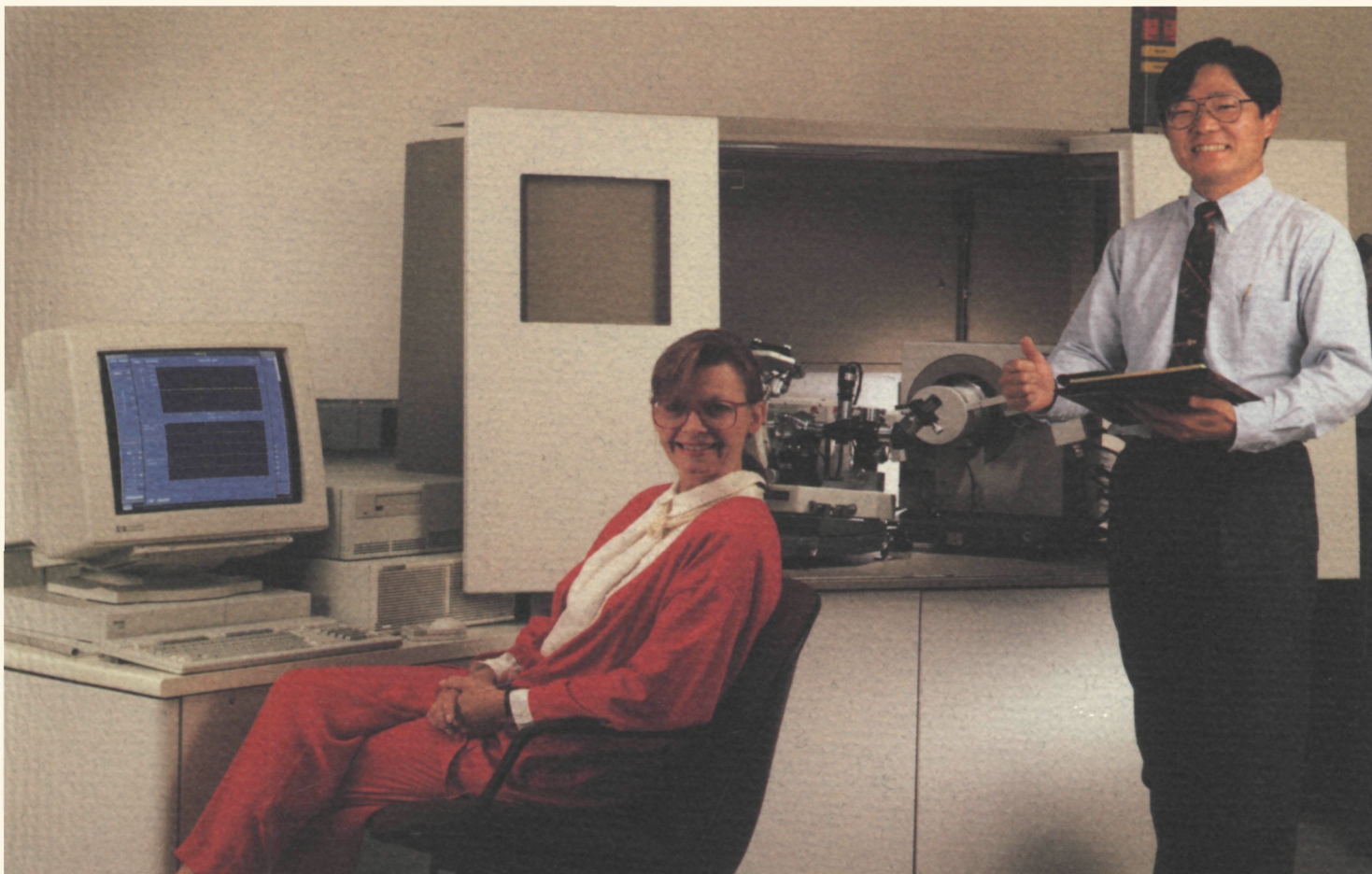

## Why are Rigaku D/MAX-1000 Diffractometers the preferred choice?

### The Reasons Are Clear.

#### Full Automation

Goniometer alignment, monochromator control, and slit selection fall into place with a single stroke of the keyboard. No other diffractometer offers such ease of use.

#### New UNIX-Based Software

Operating in a windows environment, the redesigned system boasts multiple multis. These include multi-tasking, multi-user open access, multi-windows, multi-menus and graphics screen controls, multi-benefits, and multi-cost-effectiveness.

#### Focus on Specific Needs

The D/MAX-1000 series offers a line of diffractometers that make consistent measurements unparalleled in accuracy and reproducibility. D/MAX-1000 focuses on what is needed without all the gizmos and doodads that no one ever uses

anyway. The system is designed for the precise requirements and comfort level of the individual user.

*Whatever your analysis needs, there is a D/MAX-1000 suited for you. Choose from:*

- Vertical
- Horizontal
- Theta Theta

### Contact our XRD Specialists


Learn more about the Rigaku D/MAX-1000 series and how they can provide the solutions to your analytical requirements. We develop the best diffractometers in the world, but our most important products are solutions for our customers.

**Rigaku**

#### Rigaku/USA, Inc.

Northwoods Business Park  
200 Rosewood Drive  
Danvers, MA 01923 USA  
(508) 777-2446  
FAX (508) 777-3594

# There's only one way to bring High Resolution X-ray Diffraction into your lab.

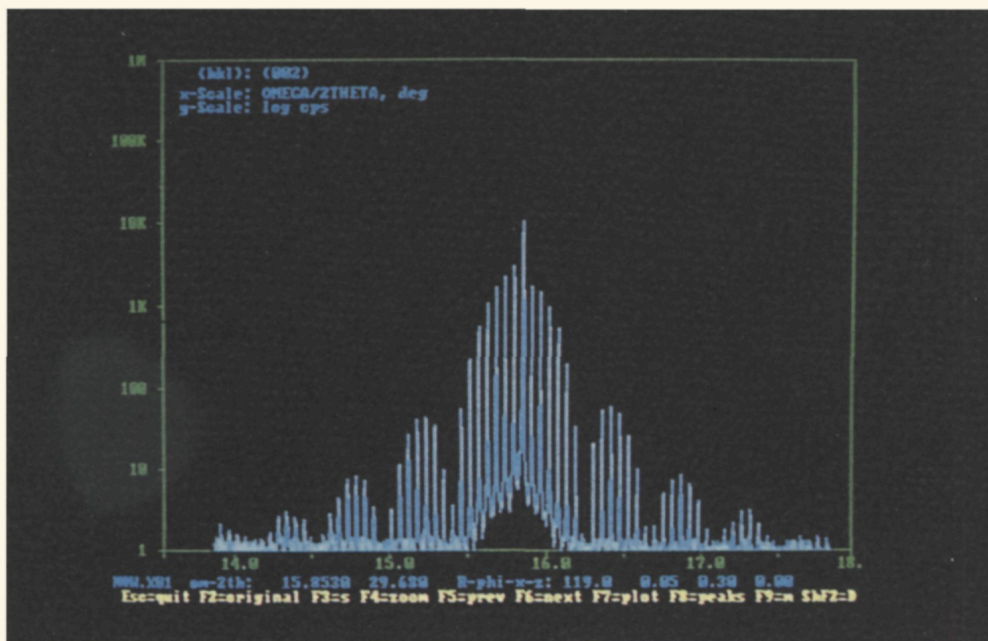


## The Philips MPD1880/HR.

This is the X-ray diffractometer that can do more for your research than any other system on the market today. It brings you that rare combination of engineering excellence in a rugged, high-performance four-crystal optical configuration that can stand up to the demands of a working lab.

The five-crystal diffraction capability provides a choice of goniometer configurations and sample cradles. That gives you more versatility in the selection of sample material to be examined or scattering geometries your protocols require. Together with the high-dynamic range sealed-xenon proportional counter, you can count on extremely sensitive measurements. And the automated sample stage allows independent theta, 2-theta, chi and phi rotation, as well as automated X-Y translation.

One of the unique features that makes the MPD1880/HR so accurate and so flexible is the Philips four-crystal monochromator. The dislocation-free, dual  $\langle 440 \rangle$  cut germanium crystals provide a wavelength definition range so small that instrumental broadening of the sample Bragg peaks is effectively removed from your analytical considerations. The user-selectable beam is the highest quality commercially available anywhere.

The choice of horizontal or vertical goniometers allows the

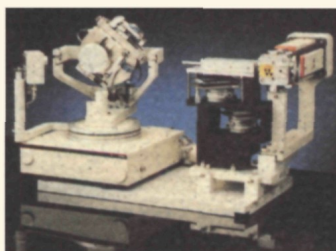

MPD1880/HR to be configured as a stand-alone system, with a dedicated X-ray tube, or as part of a multi-goniometer installation - sharing a common source with a powder, texture, or stress goniometer for complete microstructural characterization. A step size of one arc second is standard.

For fast alignments during rocking curve measurements, peak search and centering is computer controlled. A complete range of software, including

rocking curve simulation, is available so your MPD1880/HR comes with massive capability built-in.

### Features:

- Patented Philips four-crystal monochromator.
- Vertical or horizontal configurations.
- Computer-controlled sample cradles.
- X-Y stage for sample mapping.
- DOS-based software library.

For additional information contact: Philips Electronic Instruments Company, 85 McKee Drive, Mahwah, NJ 07430. Telephone (201) 529-3800.



# PHILIPS