



The Metals & Alloys Index Book for Materials Research

Materials researchers will find this new metallurgical book invaluable in their studies; its organization, selection of materials and use of alphabetical formulae is similar to that of Pearson and Villars and Calvert. The 1990 Metals & Alloys Index Book is designed to be used independently or in conjunction with the Powder Diffraction File (PDF). There are four indexes: the first two provide all materials in the Metals and Alloys PDF and the remaining two contain supporting data.

- **The Alphabetical Formula Index** brings together all entries containing a given element in alphabetical formula order. It has n entries for an n -component material and, to simplify reading the index, a straight-line format for the element being sorted. Chemical knowledge can thus be used to make identification easier.
- **The Pearson Symbol Code Index** has all entries arranged in Pearson Symbol order. Hence, the user can easily find all of the materials with a given structure, and, if desired, compare their atomic contents and lattice parameters. The latter are given for the standard Crystal Data setting to aid in structural comparisons.
- **The Common Names Index** permits cross-referencing of common metallurgical names such as austenite or cementite to the appropriate PDF data. It also has the cross-references for many mineral names that have found their way into metallurgical usage. Certain other useful names (e.g. Zintl phases) may also be found here.
- **The Strukturbericht Symbol Index** provides cross-referencing between Strukturbericht Symbols and the equivalent Pearson Symbols and Structure Prototypes. This index includes those structures for which Strukturbericht Symbols have been assigned and which are likely to be encountered in the metallurgical literature.

Price \$750.00
For information contact:

International Centre for Diffraction Data
1601 Park Lane, Swarthmore, PA 19081
USA (215) 328-9400 Telex: 847170

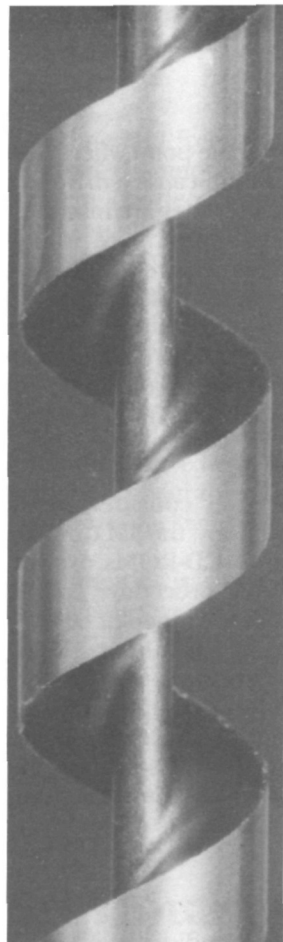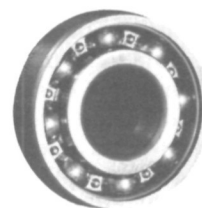

For 1992/93 from ICDD...




Set 42 PDF

In all current media
CD-ROM & Mag Tape... Microfiche ... Book
Over 2000 new &/or updated entries in
PDF-1 and PDF-2†



EDD Electron Diffraction Data Base

(NIST/Sandia/ICDD)
Crystallographic & chemical data on over
71,000 crystalline materials. †

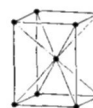


NIST Crystal Data File

Crystallographic & chemical data on over
170,000 crystalline materials*

†Available for IBM-PC compatibles, VAX & Macintosh
*Available for IBM-PC compatibles & VAX


International Centre for Diffraction Data
1601 Park Lane, Swarthmore, PA 19081
USA (215) 328-9400 Telex: 847170

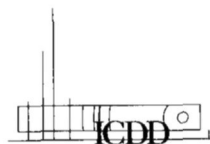

X-ray Diffraction Reference Standards and Zero-background Sample Plates Custom Designed and Built for any Application

Your first step to improved x-ray
diffraction results should be to
contact The Gem Dugout for quality
diffraction alignment standards and
zero-background plates.
And the next step is successful
x-ray diffraction results.

The Gem Dugout
1652 Princeton Drive
State College, PA 16803
(814) 865-5782

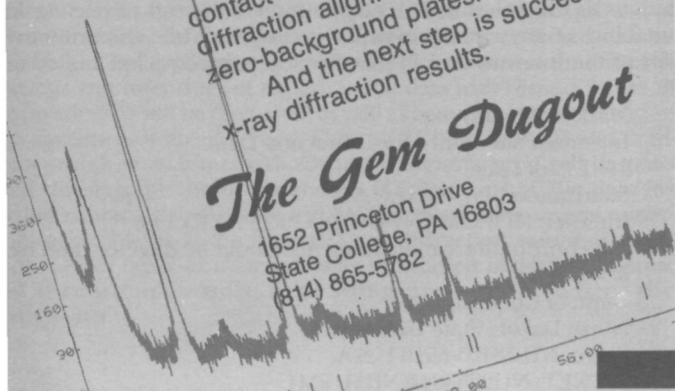

◆ ◆ ◆ ◆ ◆ ◆ ◆ ◆ ◆ ◆

Revised Sets 1-8 Databook

First in a series, this new volume was designed to replace the familiar *Black Book* version of the **Powder Diffraction File** and complement the data book version of the PDF.


Improvements in the data in this publication include:

- ◆ review through an extensive editorial process,
- ◆ evaluation by NBS* AIDS83 program,
- ◆ addition of Structure types and mineral group names,
- ◆ replacement of the ASTM phase designations by Pearson Symbol Codes,
- ◆ re-determined Quality marks
- ◆ Smith-Snyder figure of merit for each indexed entry.

On completion of the series, the user will have a uniformly formatted printed record of all data contained in the PDF.

For information and orders, contact:

International Centre for Diffraction Data,
1601 Park Lane, Swarthmore, PA 19081, U.S.A.
Phone (215) 328-9400, FAX (215) 328-2503.

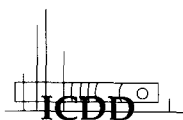

◆ ◆ ◆ ◆ ◆ ◆ ◆ ◆ ◆ ◆

ICDD Short Courses

Advance your PDF Skills

- ◆ Learn to identify/characterize crystalline substances using powder diffraction data.
- ◆ Build proficiency in interpretation of experimental powder data, using the **Powder Diffraction File**.
- ◆ Three day course includes hands-on work sessions... manual & computer search methods... comprehensive workbook.

Schedule

November 1992

Phoenix, Arizona

April 1993

Orlando, Florida


For further information contact:

Theresa Maguire

International Centre for Diffraction Data

1601 Park Lane, Swarthmore, PA 19081, U.S.A.


Phone (215) 328-9403 FAX (215) 328-2503

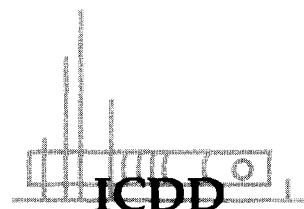

Capillary Tubes and Laboratory Supplies for X-Ray Diffraction Use

Kodak Direct Exposure Film especially developed for X-ray diffraction analysis. Available in 35mm x 20m rolls and in 50 sheet boxes.

Polaroid XR-7 High Speed X-Ray Diffraction Cassettes for Laue and Precession Photography, and Polaroid 3000 Speed Film, Type 57, in single-exposure packets.

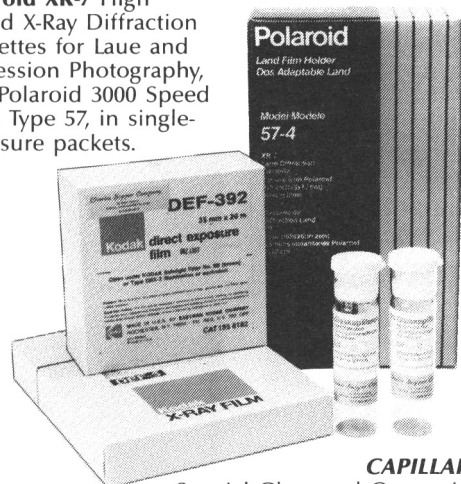

CAPILLARY TUBES

Special Glass and Quartz in 15 sizes from 0.1mm to 3.5mm O.D. Uniform 89mm length and funnel shaped at the top. 0.01mm wall thickness.

Manufacturers of Fine
X-Ray Diffraction Equip-
ment For Over 45 Years

Charles Supper Company

INCORPORATED
15 TECH CIRCLE, NATICK, MA 01760, U.S.A.
TELEPHONE: (508) 655-4610
1-800-323-9645

Consultants Directory


Powder Diffraction makes this page available for brief announcements of services offered by its readers or by other interested persons.

- Space Cost \$150 (no agency discount). This space will permit the display of a business card or other similar size text. The business card will be reproduced at no charge.
- Free space: Up to 12 lines of text immediately below the announcement.

For details write or call Mary M. Rossi,
Advertising Manager, *Powder Diffraction*. ICDD,
1601 Park Lane, Swarthmore, Pennsylvania 19081
215-328-9403 FAX 215-328-2503

Advertisers in Powder Diffraction September 1992 Volume 7 Number 3

Company	Page Number
Charles Supper Company	XIV
Charles Supper Company	XIX
Dapple Systems	XII
Fein-Marquart Associates, Inc.	I
Gem Dugout	XVII
ICDD - For 1992/1993	XVII
ICDD - Metals and Alloys	XVII
ICDD - PC-PDF CD-ROM	XV
ICDD - Revised Sets	
ICDD - Short Courses	XVIII
IC Laboratories, Inc.	XIV
inel Instrumentation Electronique	IV
Innovative Technology	X
Ital Structures	XII
KeveX Instruments	XI
Materials Data, Inc.	IX
Omni Instruments, Inc.	VI
Omni Instruments, Inc.	VII
Philips Electronic Instruments, Inc.	Back Cover
Philips Analytical	VIII
Rigaku, U.S.A.	Inside Back Cover
Scintag, Inc.	IV
Seifert X-Ray Corporation	XVI
Siemens	Inside Front Cover
Sietronics Pty. Ltd.	XIII
SPEX Industries, Inc.	II

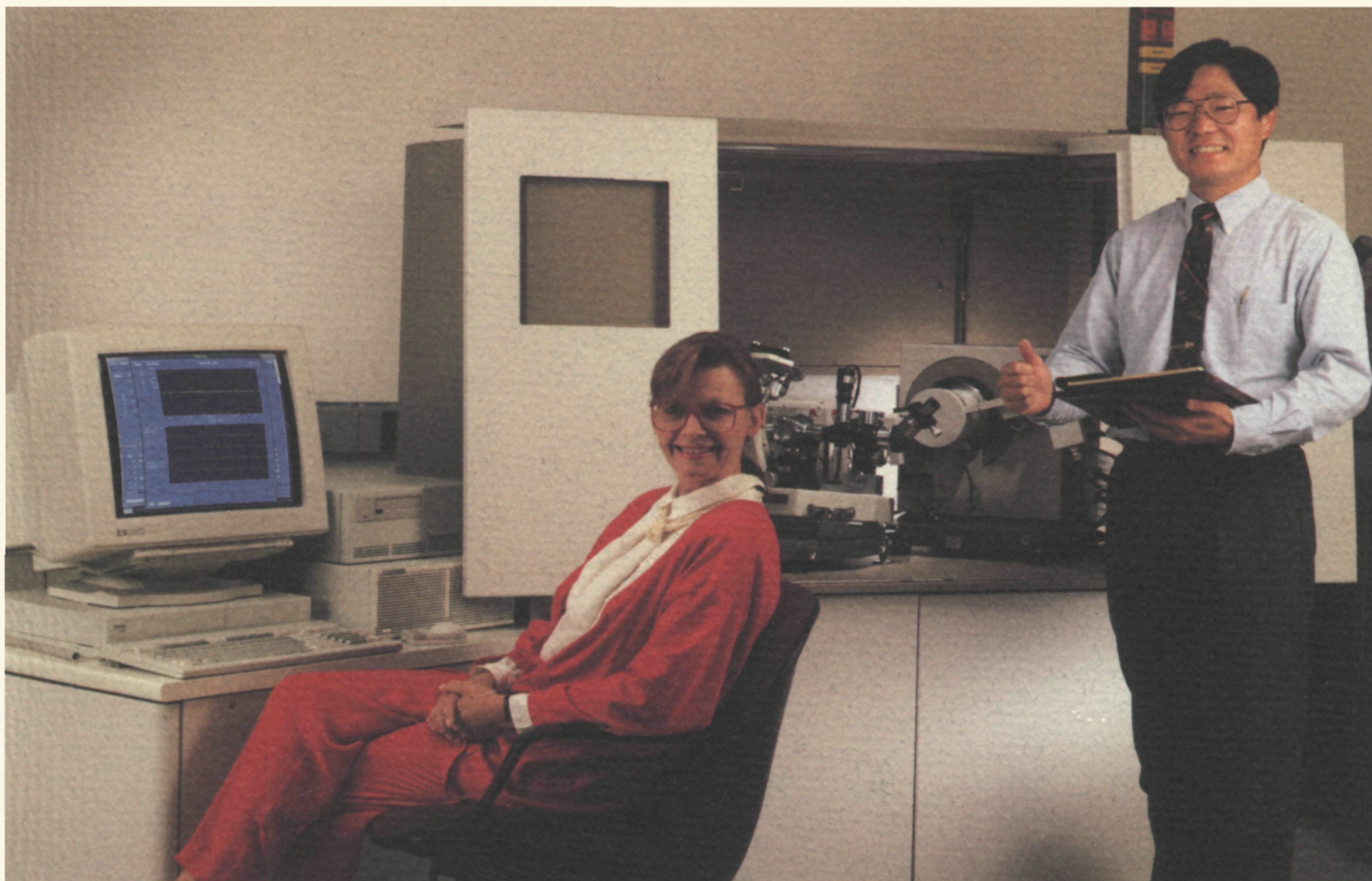

Why are Rigaku D/MAX-1000 Diffractometers the preferred choice?

The Reasons Are Clear.

Full Automation

Goniometer alignment, monochromator control, and slit selection fall into place with a single stroke of the keyboard. No other diffractometer offers such ease of use.

New UNIX-Based Software

Operating in a windows environment, the redesigned system boasts multiple multis. These include multi-tasking, multi-user open access, multi-windows, multi-menus and graphics screen controls, multi-benefits, and multi-cost-effectiveness.

Focus on Specific Needs

The D/MAX-1000 series offers a line of diffractometers that make consistent measurements unparalleled in accuracy and reproducibility. D/MAX-1000 focuses on what is needed without all the gizmos and doodads that no one ever uses

anyway. The system is designed for the precise requirements and comfort level of the individual user.

Whatever your analysis needs, there is a D/MAX-1000 suited for you. Choose from:

- Vertical
- Horizontal
- Theta Theta

Contact our XRD Specialists


Learn more about the Rigaku D/MAX-1000 series and how they can provide the solutions to your analytical requirements. We develop the best diffractometers in the world, but our most important products are solutions for our customers.

Rigaku

Rigaku/USA, Inc.

Northwoods Business Park
200 Rosewood Drive
Danvers, MA 01923 USA
(508) 777-2446
FAX (508) 777-3594

There's only one way to bring High Resolution X-ray Diffraction into your lab.

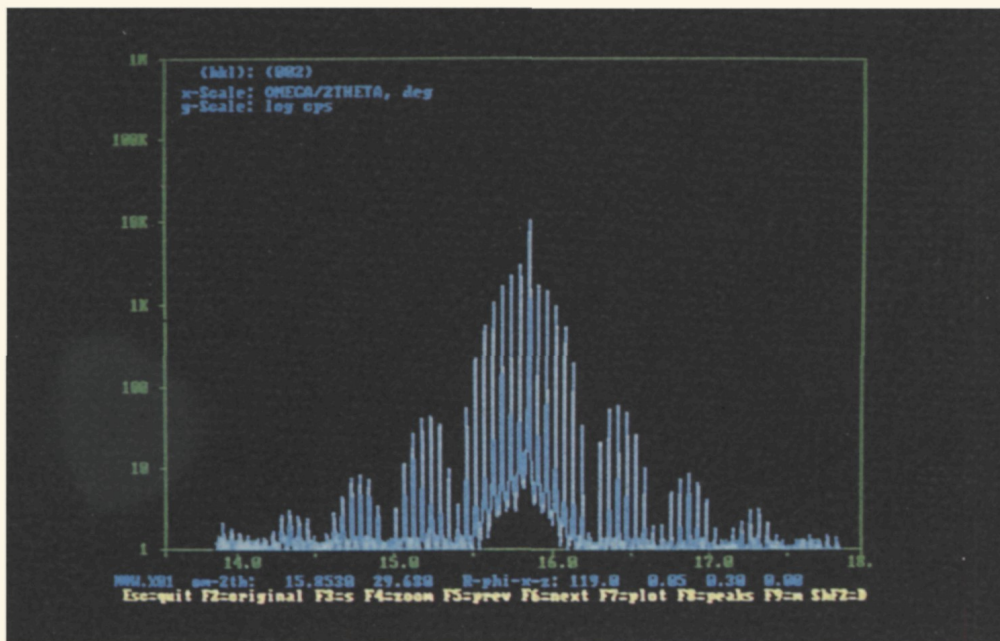


The Philips MPD1880/HR.

This is the X-ray diffractometer that can do more for your research than any other system on the market today. It brings you that rare combination of engineering excellence in a rugged, high-performance four-crystal optical configuration that can stand up to the demands of a working lab.

The five-crystal diffraction capability provides a choice of goniometer configurations and sample cradles. That gives you more versatility in the selection of sample material to be examined or scattering geometries your protocols require. Together with the high-dynamic range sealed-xenon proportional counter, you can count on extremely sensitive measurements. And the automated sample stage allows independent theta, 2-theta, chi and phi rotation, as well as automated X-Y translation.

One of the unique features that makes the MPD1880/HR so accurate and so flexible is the Philips four-crystal monochromator. The dislocation-free, dual <440> cut germanium crystals provide a wavelength definition range so small that instrumental broadening of the sample Bragg peaks is effectively removed from your analytical considerations. The user-selectable beam is the highest quality commercially available anywhere.

The choice of horizontal or vertical goniometers allows the

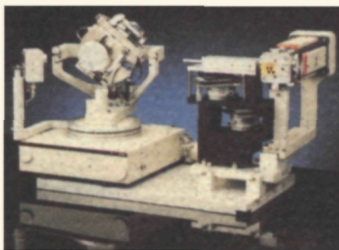

MPD1880/HR to be configured as a stand-alone system, with a dedicated X-ray tube, or as part of a multi-goniometer installation - sharing a common source with a powder, texture, or stress goniometer for complete microstructural characterization. A step size of one arc second is standard.

For fast alignments during rocking curve measurements, peak search and centering is computer controlled. A complete range of software, including rocking curve simulation, is available so your MPD1880/HR comes with massive capability built-in.

Features:

- Patented Philips four-crystal monochromator.
- Vertical or horizontal configurations.
- Computer-controlled sample cradles.
- X-Y stage for sample mapping.
- DOS-based software library.

For additional information contact: Philips Electronic Instruments Company, 85 McKee Drive, Mahwah, NJ 07430. Telephone (201) 529-3800.



PHILIPS